
UNICSKIN

YOUR EFFECTIVE SOLUTION

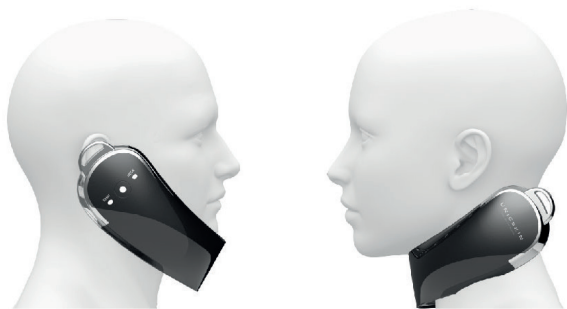UNICLED NECK & CHIN MASK

BEAUTY TECHNOLOGY

REDUCE ARRUGAS DEL CUELLO · ADELGAZA PAPADA · DEFINE ÓVALO FACIAL
LIFTS NECK · REDUCES DOUBLE CHIN · SLIMS JAWLINE

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Thank you for purchasing this product. For appropriate operation, please carefully read and keep this manual before using.

INTRODUCTION

UNICLED NECK & CHIN MASK is a diode LED light emitting device that offers a safe combination of different wavelengths. **It acts on the skin and cells, thus reducing oxidative stress and stimulating cellular energy production.** The product offers two modes: neck and chin. Depending on the mode you choose, the device will directly apply one function or another, without you having to do anything. The main benefits that you will notice will be: reduction of neck wrinkles, thinning of the double chin and jawline slimming. In addition, **the researchers argue that LED light therapy attacks adipocytes, which are cells that store fat and cause lipids to disperse. This therapy helps the body eliminate fat cells,** providing a supreme triple action when applied to the neck and chin:


- LIFTS NECK
- REDUCES DOUBLE CHIN
- SLIMS JAWLINE

USE IT 4 TIMES A WEEK. VISIBLE RESULTS IN 4 WEEKS.


WHAT TO EXPECT

Users typically see a reduction in lines and wrinkles, and a loss of fat within 2 weeks of use (4 times a week). Optimal results will start to be seen after 1 month of treatment.


LED LIGHTS WITH DEEP PENETRATION FOR BRIGHTER AND FIRMER SKIN TEXTURE:

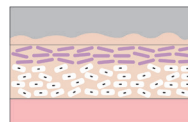Weaker elasticity.
Deeper wrinkles and fine-lines.

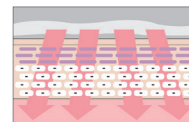RED light therapy.
Deep penetration into the skin.

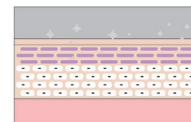Stronger elasticity.
Smoother wrinkles and fine-lines.

PACKAGE CONTENT & TECHNICAL SPECIFICATIONS

Content: device (x1), USB-C charging cable (x1) and user manual (x1).

Power: 5W

Intensity: Anti-aging 20mW/cm2 | Slimming 50mW/cm2

Voltage: 5V

Wavelength: 630nm (red) 660nm (deep red) 850nm (infrared)

Number of Lights: 112pcs LED chipset

OPERATING INSTRUCTIONS

SET UP:

- 01.** Take the **UNICLED NECK & CHIN MASK** out of the box.
- 02.** Plug the USB charging cable into the device.
- 03.** Plug USB charging cord into outlet for charging.
- 04.** Charge it for three hours (a full charge of the device takes 6 sessions).

INSTRUCTIONS FOR USE:

- 01.** Clean and dry your face and neck before starting the treatment.
- 02.** Select the desired treatment area: chin or neck. Turn on the device and let the treatment start for 20 minutes (**see more in the next section: "Benefits"**).
- 03.** The device will turn off automatically after 20 minutes.

BENEFITS

1 DEVICE, 2 MODALITIES, 3 BENEFITS. EACH MODE LASTS 20 MINUTES.

CHOOSE THE TREATMENT OF YOUR PREFERENCE:

CHIN MODE

Press the power button and choose the **"CHIN"** Mode (you will see a light on over the word). For 20 minutes, the device will act through the combination of LED lights to effectively perform 2 functions at the same time: **reduce double chin** and **jawline slimming**.

NECK MODE

Press the power button and choose **"NECK"** Mode (you will see a light on over the word). For 20 minutes, the device will act through the combination of LED lights to **effectively reduce neck wrinkles**.

IMPORTANT SAFETY INFORMATION

- Do not look directly into the lights. Placement of the device should ensure that no part of the device is covering the eyes.
- Do not store the device in direct sunlight or in a location with temperatures over 107°F/42°C as this may cause damage to the plastic coating.
- This device is designed for use on adults over 18 years of age and not recommended for use on children under the age of 12.
- Please consult your doctor before using this device if: you are or suspect that you are pregnant, nursing, if you have thyroid problems, if you are taking prescription or over-the-counter medications as these cause sensitivity to light.
- Do not place device near water or any wet surface. If device comes in contact with water, do not touch it until you first unplug it.
- Do not use device on open wounds.
- No modification of this equipment is allowed.

- Always unplug device from outlet after charging is complete.
- Do not harshly bend or abuse wire of the USB charging cord. This will void the warranty.
- Stop using device immediately if you think the device is somehow defective.

CLIENT CARE

For information regarding device usage, warranty, and additional information, please contact Customer Services by calling +34-912-797377 Monday-Friday 9:00AM-7:00PM Eastern Time. Customer Services can also be reached through e-mail at customerservices@unicskin.com.

DISCLAIMER

This product is not intended to diagnose, cure, or prevent any disease. UNICSKIN S.L. makes no claims, representations, or warranties regarding the ability of this product to cure any mental, physical, or skin condition. A qualified health professional should be consulted with regards to any condition requiring medical attention.

CLINICAL RESEARCHS

Conducted by the U.S. Department of Health & Human Services (U.S. National Library of Medicine).

In a study published in *Lasers in Surgery Medicine*, participants who received light therapy showed a significant reduction in overall circumference, including numerous parts of the body. Researchers concluded that light therapy can reduce the circumference of the specific areas of the body that are treated with natural light. **(1)**

A study conducted by Brazilian researchers in 2015 assessed the effects of light therapy and exercise on 64 obese women, ages 20 to 40. One group exercised and received red light therapy treatments, while the control group just exercised. They concluded that red light therapy and exercise together were more effective than just exercise in reducing fat mass. **(2)**

Red light therapy is safe, non-invasive, uses no chemicals or drugs, and is not associated with harmful side effects. Red light therapy simply harnesses the natural healing and rejuvenating benefits of a specific range of therapeutic natural light and delivers this targeted energy at a higher rate than the sun—without harmful UV rays. A major 2013 study on light therapy conducted by Harvard and MIT researchers praised light therapy for its “noninvasive nature and almost complete absence of side effects.” **(3)**

(1) Jackson RF, Dedo DD, Roche GC, et al. Low-level laser therapy as a non-invasive approach for body contouring: a randomized, controlled study. *Lasers in Surgery and Medicine*. Dec 2009;41(10):99-809. www.ncbi.nlm.nih.gov/pubmed/20014253

(2) Sene-Fiorese M, Duarte FO, et al. The potential of phototherapy to reduce body fat, insulin resistance and “metabolic inflexibility” related to obesity in women undergoing weight loss treatment. *Lasers Surg Med*. 2015 Oct;47(8):634-42. <https://www.ncbi.nlm.nih.gov/pubmed/26220050>

(3) Avci P, Gupta A, Sadasivam M et al. “Low-level laser (light) therapy (LLLT) in skin: stimulating, healing, restoring.” *Semin Cutan Med Surg*. 2013 March ; 32(1): 41-52. <https://www.ncbi.nlm.nih.gov/pmc/articles/PMC4126803/>

GUARANTY

Except for factory default, returns of the product are not accepted, once released and put into operation. If the packaging has not been opened, within the first 14 days, UNICSKIN will accept the return. The costs and risks of the return will be borne by the buyer. The product must be sent in the same package, along with the original proof of purchase. The shipping address is Calle José Abascal 42 2ºD, Madrid, Spain. UNICSKIN grants the legal guarantee subject to EU regulations: in case of proven factory defect, for the first 6 months, the return costs will be borne by UNICSKIN. In any case, UNICSKIN will not accept products that have suffered damage due to the customers storage and maintenance conditions.

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