# UNICSKIN

YOUR EFFECTIVE SOLUTION

## UNICLED WHITE SMILE

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## ENGLISH VERSION

Thank you for purchasing this product. For appropriate operation, please carefully read and keep this manual before using.

## INTRODUCTION

UNICLED WHITE SMILE is an innovative **Professional Mobile Teeth Whitening** to use **at-home or on-the-go. For all genders.** This kit contains everything you need to start whitening your teeth in minutes, wherever you are or whatever you are doing. Just connect the device to your Phone or Computer and use it anywhere, in the comfort of your own home, working or travelling. UNICSKIN brings you the opportunity to take advantage of the oxygenating power of the peroxide dentists use to whiten your teeth, combined with safe and effective LED Technology.

Whitening begins after your first application, and **you will see noticeably whiter, brighter teeth after 20 minutes.** The whitening gel contains ingredients for fast and effective whitening, combined with LED Accelerator Light. You can whiten your teeth in a week. The kit contains two Whitening Gel Syringes for a product use of three months. After the initial three months, you can subscribe on our website and purchase a pack of two new Whitening Gel Syringes.

## BASIC HINTS

-This product is mainly used for dental cleaning and nursing. Please carefully read the instructions before using (especially the items concerned), and then use the products correctly.

-Please correctly use it with the relevant operation instructions and requirements of this manual.

-Please immediately stop using if an exception occurs during use.

**Effectively** whiten your teeth using the most advanced teeth whitening technology and **advanced ingredients** delivering **dental grade results.** 60% of improvement in teeth whitening in the first week, 80% of improvement in the second week and 90% of improvement after 1 month of use. The same technique you'll find at the dentist, but **from the comfort and convenience of your home.** Try it out and share your results with us!

### ISSUE:

yellow teeth and other discoloration.

### TREATMENT:

ingredients on whitening gel.

### SOLUTION:


whiter and brighter teeth and smile.

**Different from other competitors, UNICSKIN's unique Whitening Gel Technology is based on an anhydrous gel, easy and safe to store for a long time; leading to higher effective whitening results.**


## CAUTION

<b>Must obey</b>	Please make sure the appearance of the product is normal before using.
<b>Must obey</b>	Please make sure to operate with the relevant requirements in this product manual.
<b>Prohibition of swallow</b>	Please rinse your mouth with warm water after using the dental cleaner and whitening gel. Don't swallow the gel.
<b>Cleaning</b>	Please clean it to avoid the accumulation of grease and other residues after using dental cleaner, keep it in clean and dry place.
	Please clean the dental cleaner with a clean wet tissue to wipe, do not rinse.
<b>Caution</b>	Please do not eat or drink alcohol within an hour of using a dental cleaner or gel pen.
	Do not look directly at the light of the dental cleaner to avoid affecting your vision.
	Please do not use if you are allergic to any of the ingredients in silicone or gel.
	If gel or vitamin E reagent contact with eyes, immediately rinse with water.
	Consult your dentist before using this product if you have dental or gum problems.
<b>Prohibition of using</b>	The following persons not suggested use: children under fourteen, pregnant or lactating women, with periodontal disease, wearing braces, dental caries or oral ulcers.
<b>Prohibition of decomposition</b>	Please do not disassemble or modify the mainframe by yourself to avoid malfunction.
<b>Prohibition of shock</b>	Please avoid falling or hitting the product to avoid malfunction.
<b>Away from children</b>	Please keep the product out of the reach of children.

## NAME OF EACH PRODUCT PART

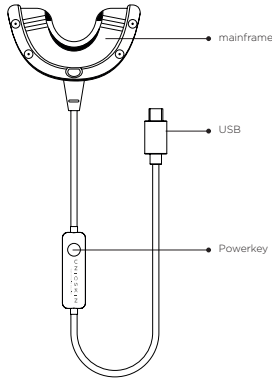
## ACCESORIES

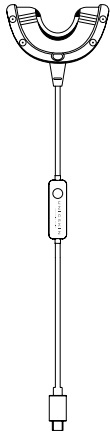
Mainframe x1


Portable bag x1

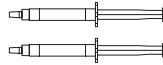Gel syringes x2

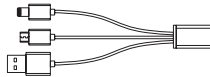Transfer interface x1

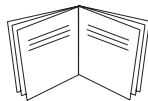
Manual x1

**Main ingredients of whitening gel:** propylene glycol, glycerin, urea, peroxide, hydrogen peroxide mixture, essential oil, triethanolamine, potassium nitrate.


## HOW TO USE

**First step:** brush your teeth normally with your regular toothpaste.

**Second step:** take out the device and clean it. Evenly apply the whitening gel on the device using the syringe.

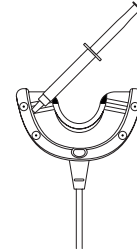**Third step:** connect the device to your smartphone or laptop. Place the device in your mouth. Press start button and leave LED blue light on.

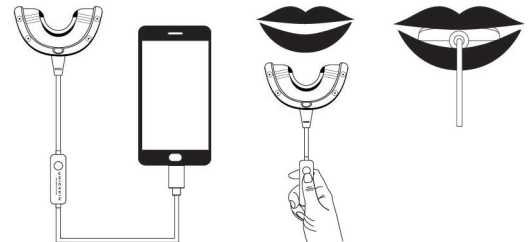**Fourth step:** after 20 minutes, the device will shut down automatically. Rinse the device and your mouth with warm water. Enjoy your brilliant teeth!

(These steps can be repeated once a week or every two weeks).

## PRODUCT STANDARD

Name of product	UNICLED WHITE SMILE
Model of product	IT02
Rated voltage	5V
Working time	Automatic shutdown after starting up for 20 minutes
Working current	MAX 240mA
Consumed power	About 1.2 W
Bare weight	About 26 g
Gel syringes	2 PCS

## WHO CAN USE AND HOW TO USE

### **Q: What kind of people is suitable to use photon dental cleaner?**

A: This product is mainly for adults with yellow teeth. This product is not recommended for children under 14 years old, pregnant or lactating women, people with periodontal disease, braces, dental caries or oral ulcers.

### **Q: How to use photon dental cleaner correctly?**

A: Use this product once a week or every two weeks.

## ADDITIONAL WAYS TO MAINTAIN YOUR WHITE SMILE

- Teeth whitening is not a permanent procedure. Maintain your white smile with regular use of the kit.
- Use an automated toothbrush.
- Practice good oral hygiene by tooth brushing, flossing and tongue cleaning.
- Regular professional dental hygiene care to maintain oral health and keep staining to a minimum.

## TRY OUR TEETH WHITENING MAGIC PEN

In addition to **UNICLED WHITE SMILE** kit, we encourage you to use our **UNICPEN WHITE SMILE**, ideal for a daily white care. Shake the pen and remove the lid. Twist the base so that the bright whitening gel flows to the bristles on the tip of the pen. Then apply the bristles evenly to your smile and focus on discoloured areas. Leave the gel on for a minute, then spit and rinse carefully.

## AFTER-SALE SERVICE

### **About internal battery**

- The internally rechargeable battery of this product, as a consumable, is beyond the guarantee range.
- The internally rechargeable battery of this product is lithium battery.

### **About the abandonment of this product**

- Please dispose this product (host/accessories) according to the local processing rules of wastes if the product is to be abandoned.

### **Consultancy center of product**

- Please consult the customer service staff if you have any questions of this product.

## GUARANTY

Except for factory default, product returns are not allowed, once they have been released and used. If the packaging has not been opened, and up to 14 days after receiving the product, UNICSKIN will accept the return. Costs and risks of the return will be borne by the buyer. The products must be sent in the same package, together with the original proof of purchase. The shipping address is Calle José Abascal, 42, 2º Right. In the event of a factory defect, UNICSKIN grants legal guarantee subject to EU regulations. During the first 6 months from the date of purchase, any return accepted by UNICSKIN will result in a replacement or a credit to the customer, at the customer's choice. In case of proven factory defect, the return costs will be borne by UNICSKIN. In any case, UNICSKIN will not accept products that have suffered damages due to the storage and maintenance conditions of the customer. The product complies with the legal and regulatory standards applicable in the European Community on the date of manufacture and distribution.

# UNICKIN

YOUR EFFECTIVE SOLUTION

UNICKIN SL. ALL RIGHTS RESERVED.

**[www.unickin.com](http://www.unickin.com)**

Contact: [customerservices@unickin.com](mailto:customerservices@unickin.com)

Calle Francisco Gervás, 14, 1ºB. 28020  
Madrid, Spain.